

## ESTRONE ELISA

EU:  IVD

CAN/USA: 

REF: CAN-E-420

Version: 9.0  
Effective: September 13, 2018

### INTENDED USE

For the direct quantitative determination of Estrone by an enzyme immunoassay in human serum.

### PRINCIPLE OF THE TEST

The principle of the following enzyme immunoassay test follows the typical competitive binding scenario. Competition occurs between an unlabelled antigen (present in standards, controls and patient samples) and an enzyme-labelled antigen (conjugate) for a limited number of antibody binding sites on the microplate. The washing and decanting procedures remove unbound materials. After the washing step, the enzyme substrate is added. The enzymatic reaction is terminated by addition of the stopping solution. The absorbance is measured on a microtiter plate reader. The intensity of the colour formed is inversely proportional to the concentration of estrone in the sample. A set of standards is used to plot a standard curve from which the amount of estrone in patient samples and controls can be directly read.

### CLINICAL APPLICATIONS

Estrone is a steroid like estradiol and estradiol, belonging to the class of estrogens. The estrogens are involved in the development of female sex organs and secondary sex characteristics. Before the ovum is fertilized the main action of the estrogens is on the growth and function of the reproductive tract in order to prepare it for the fertilized ovum.

During the follicular phase of the menstrual cycle the estrone level shows a slight increase. The production of estrone then increases markedly to peak at around day 13. The peak is of short duration and by day 16 of the cycle levels will be low. A second peak occurs at around day 21 of the cycle and if fertilization does not occur, then the production of estrone decreases.

### PROCEDURAL CAUTIONS AND WARNINGS

- Users should have a thorough understanding of this protocol for the successful use of this kit. Reliable performance will only be attained by strict and careful adherence to the instructions provided.
- Control materials or serum pools should be included in every run at a high and low level for assessing the reliability of results.
- When the use of water is specified for dilution or reconstitution, use deionized or distilled water.
- In order to reduce exposure to potentially harmful substances, gloves should be worn when handling kit reagents and human specimens.
- All kit reagents and specimens should be brought to room temperature and mixed gently but thoroughly before use. Avoid repeated freezing and thawing of reagents and specimens.

- A calibrator curve must be established for every run.
- The controls should be included in every run and fall within established confidence limits.
- Improper procedural techniques, imprecise pipetting, incomplete washing as well as improper reagent storage may be indicated when assay values for the controls do not reflect established ranges.
- When reading the microplate, the presence of bubbles in the wells will affect the optical densities (ODs). Carefully remove any bubbles before performing the reading step.
- The substrate solution (TMB) is sensitive to light and should remain colourless if properly stored. Instability or contamination may be indicated by the development of a blue colour, in which case it should not be used.
- When dispensing the substrate and stopping solution, do not use pipettes in which these liquids will come into contact with any metal parts.
- To prevent contamination of reagents, use a new disposable pipette tip for dispensing each reagent, sample, standard and control.
- Do not mix various lot numbers of kit components within a test and do not use any component beyond the expiration date printed on the label.
- Kit reagents must be regarded as hazardous waste and disposed of according to national regulations.

### LIMITATIONS

- All the reagents within the kit are calibrated for the direct determination of estrone in human serum. The kit is not calibrated for the determination of estrone in saliva, plasma or other specimens of human or animal origin.
- Do not use grossly hemolyzed, grossly lipemic, icteric or improperly stored serum.
- Any samples or control sera containing azide or thimerosal are not compatible with this kit, as they may lead to false results.
- Only calibrator A may be used to dilute any high serum samples. The use of any other reagent may lead to false results.
- The results obtained with this kit should never be used as the sole basis for a clinical diagnosis. For example, the occurrence of heterophilic antibodies in patients regularly exposed to animals or animal products has the potential of causing interferences in immunological tests. Consequently, the clinical diagnosis should include all aspects of a patient's background including the frequency of exposure to animals/products if false results are suspected.

### SAFETY CAUTIONS AND WARNINGS

#### POTENTIAL BIOHAZARDOUS MATERIAL

Human serum that may be used in the preparation of the standards and controls has been tested and found to be non-reactive for Hepatitis B surface antigen and has also been tested for the presence of antibodies to HCV and Human Immunodeficiency Virus (HIV) and found to be negative. No test method however, can offer complete assurance that HIV, HCV and Hepatitis B virus or any infectious agents are absent. The reagents should be considered a potential biohazard and handled with the same precautions as applied to any blood specimen.

#### CHEMICAL HAZARDS

Avoid contact with reagents containing TMB, hydrogen peroxide and sulfuric acid. If contacted with any of these reagents, wash with plenty of water. TMB is a suspected carcinogen.

### SPECIMEN COLLECTION AND STORAGE

Approximately 0.2 mL of serum is required per duplicate determination. Collect 4–5 mL of blood into an appropriately labelled tube and allow it to clot. Centrifuge and carefully remove the serum layer. Store at 4°C for up to 24 hours or at -10°C or lower if the analyses are to be done at a later date. Consider all human specimens as possible biohazardous materials and take appropriate precautions when handling.

### SPECIMEN PRETREATMENT

This assay is a direct system; no specimen pretreatment is necessary.

### REAGENTS AND EQUIPMENT NEEDED BUT NOT PROVIDED

- Precision pipettes to dispense 50, 100, 150 and 300 µL
- Disposable pipette tips
- Distilled or deionized water
- Plate shaker
- Microplate reader with a filter set at 450 nm and an upper OD limit of 3.0 or greater\* (see assay procedure step 10).

### REAGENTS PROVIDED

#### 1. Rabbit Anti-Estrone Antibody-Coated Break-Apart Well Microplate — Ready To Use

Contents: One 96-well (12x8) polyclonal antibody-coated microplate in a resealable pouch with desiccant.

Storage: Refrigerate at 2–8°C.

Stability: 12 months or as indicated on label.

#### 2. Estrone-Horse Radish Peroxidase (HRP) Conjugate — Ready to Use

Contents: Estrone-HRP conjugate in a protein-based buffer with a non-mercury preservative.

Volume: 14 mL/bottle

Storage: Refrigerate at 2–8°C

Stability: 12 months or as indicated on label.

#### 3. Estrone Calibrators — Ready to Use

Contents: Six vials containing estrone in a protein-based buffer with a non-mercury preservative. Prepared by spiking buffer with a defined quantity of estrone.

\* Listed below are approximate concentrations, please refer to vial labels for exact concentrations.

Calibrator	Concentration	Volume/Vial
Calibrator A	0 pg/mL	2.0 mL
Calibrator B	10 pg/mL	0.5 mL
Calibrator C	25 pg/mL	0.5 mL
Calibrator D	100 pg/mL	0.5 mL
Calibrator E	400 pg/mL	0.5 mL
Calibrator F	2000 pg/mL	0.5 mL

Storage: Refrigerate at 2–8°C

Stability: 12 months in unopened vials or as indicated on label. Once opened, the standards should be used within 14 days or aliquoted and stored frozen. Avoid multiple freezing and thawing cycles.

### 4. Controls — Ready To Use

Contents: Two vials containing estrone in a protein-based buffer with a non-mercury preservative. Prepared by spiking buffer with defined quantities of estrone. Refer to vial labels for the acceptable range.

Volume: 0.5 mL/vial

Storage: Refrigerate at 2–8°C

Stability: 12 months in unopened vial or as indicated on label. Once opened, the controls should be used within 14 days or aliquoted and stored frozen. Avoid multiple freezing and thawing cycles.

### 5. Wash Buffer Concentrate — Requires Preparation $\times 10$

Contents: One bottle containing buffer with a non-ionic detergent and a non-mercury preservative.

Volume: 50 mL/bottle

Storage: Refrigerate at 2–8°C

Stability: 12 months or as indicated on label.

Preparation: Dilute the wash buffer concentrate 1:10 in distilled or deionized water to prepare the working wash buffer. If one whole plate is to be used dilute 50 mL of the wash buffer concentrate in 450 mL of water.

### 6. TMB Substrate — Ready to Use

Contents: One bottle containing tetramethylbenzidine and hydrogen peroxide in a non-DMF or DMSO containing buffer.

Volume: 16 mL/bottle

Storage: Refrigerate at 2–8°C

Stability: 12 months or as indicated on label

### 7. Stopping Solution — Ready To Use

Contents: One bottle containing 1M sulfuric acid.

Volume: 6 mL/bottle

Storage: Refrigerate at 2–8°C

Stability: 12 months or as indicated on label.

## ASSAY PROCEDURE

Specimen Pretreatment: **None**

All reagents must reach room temperature before use. Calibrators, controls and specimen samples should be assayed in duplicate. Once the procedure has been started, all steps should be completed without interruption.

1. Prepare working solution of the wash buffer.
2. Remove the required number of well strips. Reseal the bag and return any unused strips to the refrigerator.
3. Pipette 50 µL of each calibrator, control and specimen sample into correspondingly labelled wells in duplicate.
4. Pipette 100 µL of the Estrone-HRP conjugate into each well. (We recommend using a multichannel pipette.)
5. Incubate on a plate shaker (approximately 200 rpm) for 1 hour at room temperature.
6. Wash the wells 3 times with 300 µL of diluted wash buffer per well and tap the plate firmly against absorbent paper to ensure that it is dry. (The use of a washer is recommended.)
7. Pipette 150 µL of TMB substrate into each well at timed intervals.
8. Incubate on a plate shaker for 10–15 minutes at room temperature (or until calibrator A attains dark blue colour for desired OD).
9. Pipette 50 µL of stopping solution into each well at the same timed intervals as in step 7.
10. Read the plate on a microplate reader at 450 nm within 20 minutes after addition of the stopping solution.

\* If the OD exceeds the upper limit of detection or if a 450 nm filter is unavailable, a 405 or 415 nm filter may be substituted. The optical densities will be lower, however, this will not affect the results of patient/control samples.

## CALCULATIONS


1. Calculate the mean optical density of each calibrator duplicate.
2. Draw a calibrator curve on semi-log paper with the mean optical densities on the Y-axis and the calibrator concentrations on the X-axis. If immunoassay software is being used, a 4-parameter or 5-parameter curve is recommended.
3. Calculate the mean optical density of each unknown duplicate.
4. Read the values of the unknowns directly off the calibrator curve.
5. If a sample reads more than 2000 pg/mL then dilute it with calibrator A at a dilution of no more than 1:8. The result obtained should be multiplied by the dilution factor.

## TYPICAL TABULATED DATA

Sample data only. **Do not** use to calculate results.

Calibrator	Mean OD	Value (pg/mL)
A	2.530	0
B	1.969	10
C	1.646	25
D	1.072	100
E	0.724	400
F	0.373	2000
Unknown	1.183	82.6

## TYPICAL CALIBRATOR CURVE

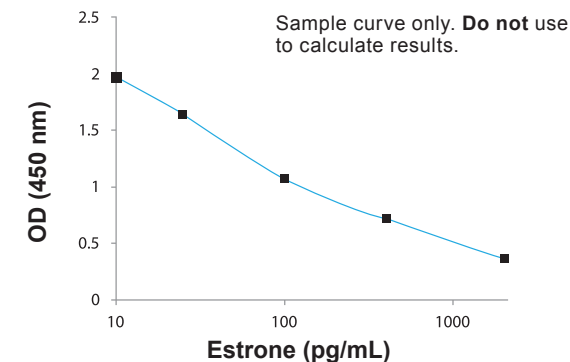

## PERFORMANCE CHARACTERISTICS

### SENSITIVITY

The limit of detection (LoD) was determined from the analysis of 60 samples of the blank and a low value sample and it was calculated as follows:

LoD =  $\mu B + 1.645\sigma B + 1.645\sigma S$ , where  $\sigma B$  and  $\sigma S$  are the standard deviation of the blank and low value sample and  $\mu B$  is the mean value of the blank.

The Limit of Detection (LoD) was determined to be **3 pg/mL**.

### SPECIFICITY (CROSS REACTIVITY)

The following compounds were tested for cross-reactivity with the Direct Estrone ELISA kit with estrone cross-reacting at 100%:

Steroid	% Cross Reactivity
Estrone	100
Estrone-3-Sulfate	4.9
17β-Estradiol	2.2
Estrone-3-Glucuronide	1.2
17β-Estradiol-3-Glucuronide	0.14

The following steroids were tested but cross-reacted at less than 0.1%: Androstenedione, Cholesterol, Corticosterone, Cortisol, Cortisone, DHEAS, Diethylstilbesterol, Estriol, 17β-Estradiol-3-Glucuronide, Estradiol-Sulfate, Progesterone, 17-OH Progesterone and Testosterone.

## INTRA-ASSAY PRECISION

Three samples were assayed twenty times each on the same calibrator curve. The results (in pg/mL) are tabulated below:

Sample	Mean	SD	CV %
1	58.1	4.76	8.2
2	239.8	20.1	8.4
3	601.3	48.7	8.1

## INTER-ASSAY PRECISION

Three samples were assayed in ten separate assays. The results (in pg/mL) are tabulated below:

Sample	Mean	SD	CV %
1	24.46	2.18	8.9
2	78.69	7.85	10.0
3	505.13	38.32	7.6

## RECOVERY

Spiked samples were prepared by adding defined amounts of estrone to three patient serum samples. The results (in pg/mL) are tabulated below:

Sample	Obs. Result	Exp. Result	Recovery %
1 Unspiked	52	-	-
+ 200	315	252	125
+ 400	557	452	120
+ 1000	1235	1052	117
2 Unspiked	75	-	-
+ 375	493	450	110
+ 750	505.23	559.81	90.3
+ 1500	712.44	794.88	89.6
3 Unspiked	720.11	-	-
+ 200	758.13	837.64	90.5
+ 400	856.46	955.17	89.7
+ 1000	1013.61	1190.24	85.1

## LINEARITY

Three patient serum samples were diluted with calibrator A. The results (in pg/mL) are tabulated below:

Sample	Obs. Result	Exp. Result	Recovery %
1	340.67	-	-
1:2	165.35	170.34	97.1
1:4	95.39	85.17	112.0
1:8	48.47	42.58	113.8
2	1086.01	-	-
1:2	508.58	543.00	93.7
1:4	232.11	271.50	85.5
1:8	114.95	135.75	84.7
3	1313.21	-	-
1:2	612.98	656.61	93.4
1:4	318.63	328.30	97.1
1:8	134.98	164.15	82.2

## COMPARATIVE STUDIES

The DBC Estrone ELISA kit (y) was compared with LC-MS/MS (x). The comparison of 19 serum samples yielded the following linear regression results:

$$y = 0.9788x - 14.656 \quad R^2 = 0.98$$

## EXPECTED NORMAL VALUES

As for all clinical assays each laboratory should collect data and establish their own range of expected normal values.

Group	Range (pg/mL)
Males	13–149
Females Pre-menopausal	21–319
Females Post-menopausal	11–95

## REFERENCES

1. Hauptmann H, et al. Concepts for the Synthesis of Biotinylated Steroids. Part 2: 17β-Estradiol Derivatives as Immunochemical Probes. *Bioconjug Chem.* 2000; 11(4):537–48.
2. Dressendörfer RA, et al. Synthesis of a Cortisol-Biotin Conjugate and Evaluation as a Tracer in an Immunoassay for Salivary Cortisol Measurement. *J Steroid Biochem Molec Biol.* 1992; 43(7):683–92.
3. Meyer HH, et al. Immunoaffinity Chromatography and a Biotin-Streptavidin Amplified Immunoassay for Sensitive and Specific Estimation of Estradiol-17β. *J Steroid Biochem.* 1990; 35(2):263–69.
4. Folan J, et al. Solid-Phase Enzymoimmunoassay of Estrone in Serum. *Clin Chem.* 1988; 34(9):1843–6.
5. Folan J, et al. Solid-Phase Enzymoimmunoassay of Estrone in Saliva. *Clin Chem.* 1989; 35(4):569–72.
6. Speroff L, et al. Hormone Biosynthesis, Metabolism, and Mechanism of Action. In: *Clinical Gynecologic Endocrinology & Infertility*, 3rd ed. Baltimore: Williams & Wilkins; 1983:1–41.
7. Kim MH, et al. Plasma Levels of Estrogen, Androgens and Progesterone During Normal and Dexamethasone-treated Cycles. *J Clin Endocrinol Metab.* 1974; 39(4):706–12.
8. Speight AC, et al. Non-Protein Bound Oestrogens in Plasma and Urinary Excretion of Unconjugated Oestrogens in Non-Pregnant Women. *J Endocrinol.* 1979; 83(3):385–91.
9. Kricka, L.J. Human Anti-Animal Antibody Interferences in Immunological Assays. *Clin Chem.* 1999; 45(7):942–56.
10. Check JH, et al. Falsely Elevated Steroidal Assay Levels Related to Heterophile Antibodies Against Various Animal Species. *Gynecol Obstet Invest.* 1995; 40(2):139–40.

**DBC-Diagnostics Biochem Canada Inc.**  
384 Neptune Crescent  
London, Ontario, Canada N6M 1A1  
Tel: (519) 681-8731  
Fax: (519) 681-8734  
e-mail: dbc@dbc-labs.com  
www.dbc-labs.com  
An ISO 13485 Registered Company

**EMERGO EUROPE**  
Prinsessegracht 20  
2514 AP The Hague  
The Netherlands

## SYMBOLS

European Conformity	In vitro diagnostic device	Consult instructions for use
Contains sufficient for <n> tests	Storage Temperature	Legal Manufacturer
Use by	Catalogue Number	Authorized representative
Lot number	Dilute 1: # Before use	